Pharma & Biotech



Efficient and reproducible results for all your transfection applications

The modular design of the Nucleofector® Platform with units for small, medium or large scale, adherent cells or cells in suspension, allows you to configure the 4D-Nucleofector® System to fit your specific research needs. Up to 3 functional units can be added to the 4D-Nucleofector® Core Unit and each is suited for different applications:

X Unit: Nucleofection® of various cell numbers in different formats

[16-well Nucleocuvette® Strips or 100 µL Nucleocuvette® Vessels]

Y Unit: Enabling adherent Nucleofection® in 24-well culture plates

LV Unit: Large-scale transfection of up to $1x10^{9}$ cells in a single run

96-well Unit: Higher throughput Nucleofection® Experiments in 96-well Nucleocuvette® Plates

There are cell type-specific Optimized Protocols or recommendations available in our knowledge center. In these Optimized Protocols the best Nucleofection® Conditions are indicated. In addition, we share our experience and knowledge for treatment of individual cell types.



Next Generation 4-D Nucleofector® Platform

Transferability of Nucleofection[®] Conditions between units and vessels

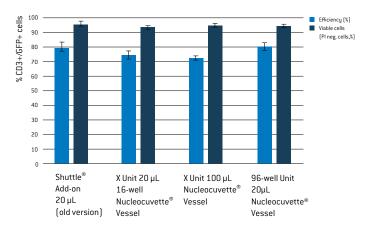


Figure 1. Transfection of human T cells with plasmid DNA (pMax-GFP® Plasmid) Fresh PBMCs were transfected with pmaxGFP® Plasmid in a 100 μ L Nucleocuvette® Vessel or a 20 μ L Nucleocuvette® Vessel using the Next Generation 4D-Nucleofector® X Unit, the 96-well Unit or the 4D-Nucleofector® Shuttle® Add-on for reference. Transfection efficiency in CD3-positive T cells was analyzed 24 hours post Nucleofection® Procedure. Plasmid data represent the mean of various independent experiments. Life cell discrimination was done with Propidium lodide.

Modular design for unmatched flexibility in scale and throughput

- All modules fully integrated in one system, including single, 16-well, 24-well and 96-well formats
- Seamless transfer of transfection conditions between different Nucleofection® Vessels
- Preconfigured system for convenient "plug & play" setup in less than two minutes
- Easy setup and management of experiments via intuitive software and touchscreen display
- High transfection efficiency, up to 90% for plasmid DNA and 99% for oligonucleotides
- Provides excellent preservation of physiological status and viability of transfected cells
- Optimized transfection protocols for >750 cell types

Ordering information

oracing morniage.		
Cat. no.	Product name	Capacity
AAF-1003B	4D-Nucleofector® Core Unit	1 Unit
AAF-1003X	4D-Nucleofector® X Unit	1 Unit
AAF-1003Y	4D-Nucleofector® Y Unit	1 Unit
AAF-1003S	4D-Nucleofector® 96-well Unit	1 Unit
AWA-3001-B	4D-Nucleofector® Core Unit Guarantee Extension	For 1 year
AWA-3001-X	4D-Nucleofector® X Unit Guarantee Extension	For 1 year
AWA-3001-Y	4D-Nucleofector® Y Unit Guarantee Extension	For 1 year
AWA-3001-96	4D-Nucleofector® 96-well Unit Guarantee Extension	For 1 uear

Contact us

North America

Customer Service: +1 800 638 8174 (toll free) Fax: +1 301 845 8338

order.us@lonza.com

Scientific Support: +1 800 521 0390 (toll free)

scientific.support@lonza.com

Europe

Customer Service: +32 87 321 611

order.europe@lonza.com

Scientific Support: +32 87 321 611 scientific.support.eu@lonza.com

International

Lonza Cologne GmbH – 50829 Cologne Germany

For research use only. Not for use in diagnostic procedures. The Nucleofector® Technology is covered by patent and/or patent pending rights owned by the Lonza Group Ltd or its affiliates. All trademarks belong to Lonza, registered in USA, EU or CH or to third party owners and used only for informational purposes. All third party copyrights have been reproduced with permission or belong to Lonza. The information contained herein is believed to be correct and corresponds to the latest state of scientific and technical knowledge. However, no warranty is made, either expressed or implied, regarding its accuracy or the results to be obtained from the use of such information and no warranty is expressed or implied concerning the use of these products. The buyer assumes all risks of use and/or handling. Any user must make his own determination and satisfy himself that the products supplied by Lonza Group Ltd or its affiliates and the information and recommendations given by Lonza Group Ltd or its affiliates are (i) suitable for intended process or purpose, (ii) in compliance with environmental, health and safety regulations, and (iii) will not infringe any third party's intellectual property rights. The user bears the sole responsibility for determining the existence of any such third party rights, as well as obtaining any necessary licenses. For more details: www.lonza.com/legal.

© 2021 Lonza. All rights reserved. CD-TS117 08/21